

More Chemistries, More Choices for Your Toughest Separation Challenges

Agilent InfinityLab Poroshell 120 columns



Efficiently separate the widest variety of columns

Based on superficially porous particle technology, InfinityLab Poroshell 120 columns feature a solid silica core and a porous outer layer. Compared to traditional totally porous particles of the same (or similar) size, Poroshell particles deliver higher chromatographic efficiencies and enable fast, high-resolution separations.

The InfinityLab Poroshell 120 family includes 20 chemistries. This ordering guide will help you choose the right particle size and bonded phase to achieve optimal separation for your samples.

Which particle is best for my method?

InfinityLab Poroshell 120 columns are available in three particle dimensions. Choose the size that best fits your separation needs and is most compatible with your LC systems.

UHPLC		1.9 µm: highest UHPLC performance – Maximum pressure: 1300 bar – Ideal for: Agilent 1290 Infinity II LC
HPLC or UHPLC		2.7 µm: UHPLC performance at lower pressures – Maximum pressure: 600 bar (unless otherwise noted) – Ideal for: Agilent 1260 Infinity II LC or Agilent 1260 Infinity II Prime LC
HPLC		4 µm: improved HPLC performance – Maximum pressure: 600 bar – Ideal for: Agilent 1220 Infinity II LC

What column ID and length should I choose?

Format	Comment
Column ID	4.6 mm for legacy methods 3.0 mm for lower solvent use than 4.6 mm 2.1 mm for lowest solvent use and MS applications
Column length	Shorter 30 to 100 mm for fastest separations Longer 150 to 250 mm for increased resolution

InfinityLab Poroshell 120 column specifications

InfinityLab Poroshell 120	Pore Size	Temperature Limit	pH Range	Endcapped	Carbon Load	Surface Area	USP Designation
EC-C18	120 Å	60 °C	2.0–8.0	Yes	10%	130 m ² /g	L1
EC-C8	120 Å	60 °C	2.0–8.0	Yes	5%	130 m ² /g	L7
Aq-C18	120 Å	90 °C	1.0–8.0	Yes	Proprietary	130 m ² /g	L1
SB-C18	120 Å	90 °C	1.0–8.0	No	9%	130 m ² /g	L1
SB-C8	120 Å	80 °C	1.0–8.0	No	5.5%	130 m ² /g	L7
CS-C18	100 Å	90 °C	1.0–11.0	Yes	Proprietary	95 m ² /g	L1
HPH-C18	100 Å	60 °C	2.0–11.0	Yes	Proprietary	95 m ² /g	L1
HPH-C8	100 Å	60 °C	2.0–11.0	Yes	Proprietary	95 m ² /g	L7
Bonus-RP	120 Å	60 °C	2.0–8.0	Yes	9.5%	130 m ² /g	L60
PFP	120 Å	60 °C	2.0–8.0	Yes	5.1%	130 m ² /g	L43
Phenyl-Hexyl	120 Å	60 °C	2.0–8.0	Yes	9%	130 m ² /g	L11
SB-Aq	120 Å	80 °C	1.0–8.0	No	Proprietary	130 m ² /g	L96
EC-CN	120 Å	60 °C	2.0–8.0	Yes	3.5%	130 m ² /g	L10
HILIC-Z	100 Å	80 °C	2.0–12.0	No	Proprietary	95 m ² /g	L114
HILIC	120 Å	60 °C	1.0–8.0	No	NA	130 m ² /g	L3
HILIC-OH5	120 Å	45 °C	1.0–7.0	Proprietary	Proprietary	130 m ² /g	L86
Chiral-V	120 Å	45 °C	2.5–7.0	Proprietary	Proprietary	130 m ² /g	L88
Chiral-T	120 Å	45 °C	2.5–7.0	Proprietary	Proprietary	130 m ² /g	L63
Chiral-CD	120 Å	45 °C	3.0–7.0	Proprietary	Proprietary	130 m ² /g	L45
Chiral-CF	120 Å	45 °C	3.0–7.0	Proprietary	Proprietary	130 m ² /g	NA

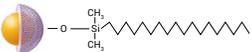
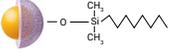
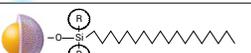
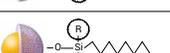
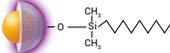
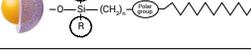
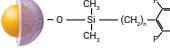
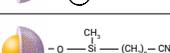
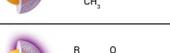
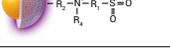
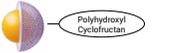
What if my methods were developed on fully porous columns?

InfinityLab Poroshell chemistries are aligned with traditional ZORBAX chemistries—making it easy to transfer your methods from fully porous to superficially porous particle columns.

InfinityLab Poroshell Chemistry	Aligned Chemistry
InfinityLab Poroshell 120 EC-C18	ZORBAX Eclipse Plus C18
InfinityLab Poroshell 120 EC-C8	ZORBAX Eclipse Plus EC-C8
InfinityLab Poroshell 120 Phenyl-Hexyl	ZORBAX Eclipse Plus Phenyl-Hexyl
InfinityLab Poroshell 120 SB-C18	ZORBAX StableBond SB-C18
InfinityLab Poroshell 120 SB-C8	ZORBAX StableBond SB-C8
InfinityLab Poroshell 120 Bonus-RP	ZORBAX Bonus-RP
InfinityLab Poroshell 120 SB-Aq	ZORBAX StableBond SB-Aq
InfinityLab Poroshell 120 EC-CN	ZORBAX Eclipse XDB-CN
InfinityLab Poroshell 120 HILIC	ZORBAX HILIC Plus

Which bonded phase will help me perfect my separations?

With 20 different chemistries, InfinityLab Poroshell 120 columns provide a range of selectivities that makes your method development fast and easy.

InfinityLab Poroshell 120	Chemistry	Particle Sizes	Benefits and Applications
EC-C18		1.9, 2.7, 4 μm	General purpose Excellent peak shape and efficiency for acids, bases, and neutrals
EC-C8		1.9, 2.7, 4 μm	General purpose Lower retention of hydrophobic analytes vs. C18
Aq-C18		2.7 μm	Enhanced retention for challenging polar compounds while also separating non-polar analytes. 100% aqueous mobile phase compatibility and low pH stability
SB-C18		1.9, 2.7, 4 μm	Low pH Excellent stability and peak shape in highly acidic conditions
SB-C8		2.7 μm	Low pH Excellent stability at low pH Lower retention of hydrophobic analytes vs. C18
CS-C18		2.7 μm	Alternate selectivity Improved peak shape and sample capacity for basic compounds with low ionic strength mobile phases High pH capable
HPH-C18		1.9, 2.7, 4 μm	High pH capable Robust performance and long lifetimes Improved retention, resolution, and peak shape of basic compounds
HPH-C8		2.7, 4 μm	High pH capable Robust performance and long lifetimes Lower retention of hydrophobic analytes vs. C18
Bonus-RP		2.7 μm	Alternative selectivity to C18 Improved peak shape for basic compounds, stable in 100% aqueous conditions
PFP		1.9, 2.7, 4 μm	Alternative selectivity Excellent peak shape for polar and nonpolar analytes Unique selectivity for aromatic and halogenated compounds
Phenyl-Hexyl		1.9, 2.7, 4 μm	Alternative selectivity with aromatic groups Highly nonpolar bonded phase takes advantage of pi-pi interactions
SB-Aq		1.9, 2.7, 4 μm	Alternative selectivity Excellent peak shape and retention of polar compounds using reversed-phase LC Exceptional stability under high-aqueous conditions, including 100% water
EC-CN		2.7 μm	Alternative selectivity Use in reversed phase for alternative selectivity of polar and midpolar compounds Use in normal phase for excellent peak shape and retention of nonpolar analytes
HILIC-Z		1.9, 2.7, 4 μm	Polar analytes Excellent retention of highly polar or charged compounds by HILIC Rugged performance at high pH or high temperature
HILIC		1.9, 2.7, 4 μm	Polar analytes Excellent retention of polar compounds by HILIC
HILIC-OH5		2.7 μm	Polar analytes Fructan bonded phase offers alternative selectivity to other HILIC phases
Chiral-V		2.7 μm	Chiral separations Amines, profens, and complex basic and neutral compounds Reversed phase, polar ionic, normal phase, or polar organic modes
Chiral-T		2.7 μm	Chiral separations Beta blockers, hydroxyl acids, amino acids, profens, benzodiazepines, and hydantoins Reversed phase, polar ionic, or polar organic modes
Chiral-CD		2.7 μm	Chiral separations Stimulants, fungicides, and protected amino acids Reversed phase or polar organic modes
Chiral-CF		2.7 μm	Chiral separations Primary amines Polar organic or normal phase modes

HILIC-OH5 and all four chiral phases have a pressure limit of 400 bar.

Ordering information

InfinityLab Poroshell 120 1.9 µm columns: Highest UHPLC performance

Size (mm)	EC-C18	EC-C8	Phenyl-Hexyl	HPH-C18	PFP	HILIC	HILIC-Z	SB-C18	SB-Aq
3.0 x 150	693675-302	693675-306	693675-312	693675-502	693675-308	693675-301		683675-302	
3.0 x 100	695675-302	695675-306	695675-312	695675-502	695675-308	695675-301		685675-302	
3.0 x 50	699675-302	699675-306	699675-312	699675-502	699675-308	699675-301		689675-302	
3.0 X 30	691775-302			691775-502					
2.1 x 150	693675-902	693675-906	693675-912	693675-702	693675-408	693675-901	683675-924	683675-902	683675-914
2.1 x 100	695675-902	695675-906	695675-912	695675-702	695675-408	695675-901	685675-924	685675-902	685675-914
2.1 x 50	699675-902	699675-906	699675-912	699675-702	699675-408	699675-901	689675-924	689675-902	689675-914

InfinityLab Poroshell 120 1.9 µm columns have a 1300 bar/18,850 psi pressure limit.

Except for 1220 InfinityLab II, all new instruments can be equipped with column ID Tag Readers to monitor tagged column usage status for ease of use, traceability, and safety. Currently, selected InfinityLab Poroshell 120 columns come pre-tagged. To order columns with a column ID tag, purchase part number with "T" at the end (example 690970-902T). If there is not a pre-tagged column option, please reach out to your CSD sales representative for further assistance.

If you need any support on applications, please contact CSD LC column support in your region.

InfinityLab Poroshell 120 1.9 µm UHPLC Fast Guards: Extend column lifetime even further

Size (mm)	EC-C18	EC-C8	Phenyl-Hexyl	HPH-C18	PFP	HILIC	SB-C18
3.0 x 5	823750-940	823750-941	823750-943	823750-945	823750-942	823750-944	823750-950
2.1 x 5	821725-940	821725-941	821725-943	821725-945	821725-942	821725-944	821725-948

Guards supplied as 3/pkg.

InfinityLab Poroshell 120 2.7 µm chiral columns: High-resolution separations of chiral compounds

Size (mm)	Chiral-CF	Chiral-CD	Chiral-V	Chiral-T
2.1 x 50	689775-609	689775-607	689775-604	689775-603
2.1 x 100	685775-609	685775-607	685775-604	685775-603
2.1 x 150	683775-609	683775-607	683775-604	683775-603
4.6 x 50	689975-609	689975-607	689975-604	689975-603
4.6 x 100	685975-609	685975-607	685975-604	685975-603
4.6 x 150	683975-609	683975-607	683975-604	683975-603

Rated to 400 bar. InfinityLab Poroshell 120 chiral columns are not available with preprogrammed column ID tags.

InfinityLab Poroshell 120 2.7 µm HILIC columns: Best for polar separations

Size (mm)	HILIC	HILIC-OH5	HILIC-Z	HILIC-Z (PEEK-lined)
2.1 x 50	699775-901	689775-601	689775-924	679775-924
2.1 x 100	695775-901	685775-601	685775-924	675775-924
2.1 x 150	693775-901	683775-601	683775-924	673775-924
3.0 x 50	699975-301		689975-324	
3.0 x 100	695975-301		685975-324	
3.0 x 150	693975-301		683975-324	
4.6 x 50	699975-901	689975-601	689975-924	
4.6 x 100	695975-901	685975-601	685975-924	
4.6 x 150	693975-901	683975-601	683975-924	

InfinityLab Poroshell 120 HILIC-OH5 has a 400 bar maximum pressure limit and is not available with a preprogrammed column ID.



For more information and helpful tools for chiral testing, including the chiral application e-notebook with more than 200 application examples, visit www.agilent.com/chem/poroshell-120-chiral

InfinityLab Poroshell 120 2.7 µm columns: UHPLC performance at lower pressures

Size (mm)	EC-C18	EC-C8	Aq-C18	SB-C18	SB-C8	CS-C18	CS-C18 PEEK-lined	HPH-C18	HPH-C8
4.6 x 150	693975-902	693975-906	693975-742	683975-902	683975-906	693975-942		693975-702	693975-706
4.6 x 100	695975-902	695975-906	695975-742	685975-902	685975-906	695975-942		695975-702	695975-706
4.6 x 75	697975-902	697975-906		687975-902					
4.6 x 50	699975-902	699975-906	699975-742	689975-902	689975-906	699975-942		699975-702	699975-706
4.6 x 30	691975-902	691975-906		681975-902					
3.0 x 150	693975-302	693975-306		683975-302	683975-306	693975-342		693975-502	693975-506
3.0 x 150	693575-302*		693675-742*						
3.0 x 100	695975-302	695975-306		685975-302	685975-306	695975-342		695975-502	695975-506
3.0 x 100	695575-302*		695675-742*						
3.0 x 75	697975-302	697975-306		687975-302					
3.0 x 50	699975-302	699975-306	699675-742	689975-302	689975-306	699975-342		699975-502	699975-506
3.0 x 30	691975-302	691975-306		681975-302					
2.1 x 150	693775-902	693775-906		683775-902	683775-906	693775-942		693775-702	693775-706
2.1 x 150	693575-902*		693775-742*						
2.1 x 100	695775-902	695775-906		685775-902	685775-906	695775-942	675775-942	695775-702	695775-706
2.1 x 100	695575-902*		695775-742*						
2.1 x 75	697775-902	697775-906		687775-902					
2.1 x 50	699775-902	699775-906	699775-742	689775-902	689775-906	699775-942	679775-942	699775-702	699775-706
2.1 x 30	691775-902	691775-906	691775-742	681775-902					
Size (mm)	Phenyl-Hexyl	Bonus-RP	PFP	SB-Aq	EC-CN				
4.6 x 150	693975-912	693968-901	693975-408	683975-914	693975-905				
4.6 x 100	695975-912	695968-901	695975-408	685975-914	695975-905				
4.6 x 50	699975-912	699968-901	699975-408	689975-914	699975-905				
3.0 x 150	693975-312	693968-301	693975-308	683975-314	693975-305				
3.0 x 100	695975-312	695968-301	695975-308	685975-314	695975-305				
3.0 x 50	699975-312	699968-301	699975-308	689975-314	699975-305				
2.1 x 150	693775-912	693768-901	693775-408	683775-914	693775-905				
2.1 x 100	695775-912	695768-901	695775-408	685775-914	695775-905				
2.1 x 50	699775-912	699768-901	699775-408	689775-914	699775-905				

InfinityLab Poroshell 120 2.7 µm columns have a 600 bar/9,000 psi pressure limit unless otherwise noted.

*These columns have a maximum pressure limit of 1,000 bar/14,500 psi.

Except for 1220 InfinityLab II, all new instruments can be equipped with column ID Tag Readers to monitor tagged column usage status for ease of use, traceability, and safety. Currently, selected InfinityLab Poroshell 120 columns come pre-tagged. To order columns with a column ID tag, purchase part number with "T" at the end (example 690970-902T). If there is not a pre-tagged column option, please reach out to your CSD sales representative for further assistance.

If you need any support on applications, please contact CSD LC column support in your region.

InfinityLab Poroshell 120 2.7 µm Fast Guards: Extend column lifetime even further

Size (mm)	EC-C18	EC-C8	Aq-C18	SB-C18	SB-C18	CS-C18	HPH-C18	
4.6 x 5	820750-911	820750-913	820750-942	820750-912	820750-923	820750-939	820750-928	
3.0 x 5	823750-911	823750-913	823750-953	823750-912	823750-923	823750-949	821725-928	
2.1 x 5	821725-911	821725-913	821725-955	821725-912	821725-923	821725-953	823750-928	
Size (mm)	HPH-C8	Phenyl-Hexyl	Bonus-RP	PFP	SB-Aq	EC-CN	HILIC	HILIC-Z
4.6 x 5	820750-922	820750-914	820750-925	820750-915	820750-924	820750-927	820750-926	820750-933
3.0 x 5	823750-922	823750-914	823750-925	823750-915	823750-924	823750-927	823750-926	823750-948
2.1 x 5	821725-922	821725-914	821725-925	821725-915	821725-924	821725-927	821725-926	821725-947

Guards supplied as 3/pk.

InfinityLab Poroshell 120 4 µm columns: Improved HPLC performance

Size (mm)	EC-C18	EC-C8	Phenyl-Hexyl	HPH-C18	HPH-C8	PFP	HILIC	HILIC-Z	SB-C18	SB-Aq
4.6 x 250	690970-902	690970-906	690970-912	690970-702	690970-706	690970-408	690970-901	680970-924	680970-902	680970-914
4.6 x 150	693970-902	693970-906	693970-912	693970-702	693970-706	693970-408	693970-901	683970-924	683970-902	683970-914
4.6 x 100	695970-902	695970-906	695970-912	695970-702	695970-706	695970-408	695970-901	685970-924	685970-902	685970-914
4.6 x 50	699970-902	699970-906	699970-912	699970-702	699970-706	699970-408	699970-901	689970-924	689970-902	689970-914
3.0 x 250	690970-302	690970-306	690970-312	690970-502	690970-506	690970-308	690970-301		680970-302	
3.0 x 150	693970-302	693970-306	693970-312	693970-502	693970-506	693970-308	693970-301		683970-302	
3.0 x 100	695970-302	695970-306	695970-312	695970-502	695970-506	695970-308	695970-301		685970-302	
3.0 x 50	699970-302	699970-306	699970-312	699970-502	699970-506	699970-308	699970-301		689970-302	
2.1 x 250	650750-902	650750-906	650750-912	690770-702	690770-706	650750-408	650750-901		640750-902	
2.1 x 150	693770-902	693770-906	693770-912	693770-702	693770-706	693770-408	693770-901		683770-902	
2.1 x 100	695770-902	695770-906	695770-912	695770-702	695770-706	695770-408	695770-901		685770-902	
2.1 x 50	699770-902	699770-906	699770-912	699770-702	699770-706	699770-408	699770-901		689770-902	

InfinityLab Poroshell 120 4 µm columns have a 600 bar/9,000 psi pressure limit.

Except for 1220 InfinityLab II, all new instruments can be equipped with column ID Tag Readers to monitor tagged column usage status for ease of use, traceability, and safety. Currently, selected InfinityLab Poroshell 120 columns come pre-tagged. To order columns with a column ID tag, purchase part number with "T" at the end (example 690970-902T). If there is not a pre-tagged column option, please reach out to your CSD sales representative for further assistance.

If you need any support on applications, please contact CSD LC column support in your region.

InfinityLab Poroshell 120 4 µm Fast Guards: Extend column lifetime even further

Size (mm)	EC-C18	HPH-C18	HPH-C8	SB-C18
4.6 x 5	820750-916	820750-928	820750-929	820750-934
3.0 x 5	823750-916	823750-928	823750-929	823750-951
2.1 x 5	821725-916	821725-928	821725-929	821725-949

Guards supplied as 3/pk.

InfinityLab Poroshell 120 preparative 4 µm columns

Size (mm)	SB-C18	HPH-C18
21.2 x 50	670050-902	670050-702
21.2 x 150	670150-902	670150-702



Packed with 4 µm superficially porous particles (SPP), Agilent InfinityLab Poroshell 120 preparative LC columns provide the high-resolution separations you need to purify your complex samples. They increase your throughput by maintaining this performance at higher flow rates, saving you time from method development to purification. So, you can successfully complete your purification goals faster and increase lab productivity.

www.agilent.com/chem/prepcolumns

Method validation kits

Kits contain three HPLC columns with the same column type (bonded phase, particle size, configuration), but from different manufacturing lots.

Size (mm)	Particle Size (µm)	EC-C18	EC-C8	Phenyl-Hexyl	SB-C18	SB-C8	SB-Aq	Bonus-RP
4.6 x 150	2.7	693975-902K	693975-906K	693975-912K	683975-902K	683975-906K	683975-914K	693968-901K
4.6 x 100	2.7	695975-902K	695975-906K	695975-912K	685975-902K	685975-906K	685975-914K	695968-901K
4.6 x 50	2.7	699975-902K	699975-906K	699975-912K	689975-902K	689975-906K	689975-914K	699968-901K
3.0 x 150	2.7	693975-302K	693975-306K	693975-312K	683975-302K	683975-306K	683975-314K	693968-301K
3.0 x 100	2.7	695975-302K	695975-306K	695975-312K	685975-302K	685975-306K	685975-314K	695968-301K
3.0 x 50	2.7	699975-302K	699975-306K	699975-312K	689975-302K	689975-306K	689975-314K	699968-301K
2.1 x 150	2.7	693775-902K	693775-906K	693775-912K	683775-902K	683775-906K	683775-914K	693768-901K
2.1 x 100	2.7	695775-902K	695775-906K	695775-912K	685775-902K	685775-906K	685775-914K	695768-901K
2.1 x 50	2.7	699775-902K	699775-906K	699775-912K	689775-902K	689775-906K	689775-914K	699768-901K
Size (mm)	Particle Size (µm)	EC-C18	EC-C8	Phenyl-Hexyl	PPF	HILIC		
4.6 x 250	4	690970-902K	690970-906K	690970-912K	690970-408K	690970-901K		
4.6 x 150	4	693970-902K	693970-906K	693970-912K	693970-408K	693970-901K		
4.6 x 100	4	695970-902K	695970-906K	695970-912K	695970-408K	695970-901K		
4.6 x 50	4	699970-902K	699970-906K	699970-912K	699970-408K	699970-901K		
3.0 x 250	4	690970-302K	690970-306K	690970-312K	690970-308K	690970-301K		
3.0 x 150	4	693970-302K	693970-306K	693970-312K	693970-308K	693970-301K		
3.0 x 100	4	695970-302K	695970-306K	695970-312K	695970-308K	695970-301K		
3.0 x 50	4	699970-302K	699970-306K	699970-312K	699970-308K	699970-301K		
2.1 x 250	4	650750-902K	650750-906K	650750-912K	650750-408K	650750-901K		
2.1 x 150	4	693770-902K	693770-906K	693770-912K	693770-408K	693770-901K		
2.1 x 100	4	695770-902K	695770-906K	695770-912K	695770-408K	695770-901K		
2.1 x 50	4	699770-902K	699770-906K	699770-912K	699770-408K	699770-901K		

Don't see a validation kit for your column part number?

Validation kits can be assembled for most columns in the Agilent LC column portfolio. Simply contact your Agilent sales specialist or send your request to the Agilent customer center in your country, indicating the part number for which you need a validation kit. An Agilent representative will send you a quote within one to two business days. www.agilent.com/chem/contactus

Learn more:

www.agilent.com/chem/poroshell-120

U.S. and Canada

1-800-227-9770

agilent_inquiries@agilent.com

Europe

info_agilent@agilent.com

Asia Pacific

inquiry_lsca@agilent.com

www.agilent.com/chem/small-molecule-columns

DE41930572

This information is subject to change without notice.

© Agilent Technologies, Inc. 2022
Printed in the USA, December 22, 2022
5991-9123EN

